

**LINK**<sup>®</sup> 

*LinkSymphoKnee*  
Built on TRUST 




*“The LinkSymphoKnee PS knee was conceived to achieve predictable femoral rollback. The components also allow the use of a PS+ insert to seamlessly address intraoperative ligament instability, without a detrimental effect on range of motion. Our early outcomes have been very favorable with excellent postoperative motion, function and pain relief.”*

**– Alejandro G. Della Valle MD Hospital for Special Surgery, New York City, NY<sup>1</sup>**

<sup>1</sup>Dr. Della Valle performed the first United States implantation of the LSK at HSS in August 2021.



*LinkSymphoKnee*  
Built on TRUST

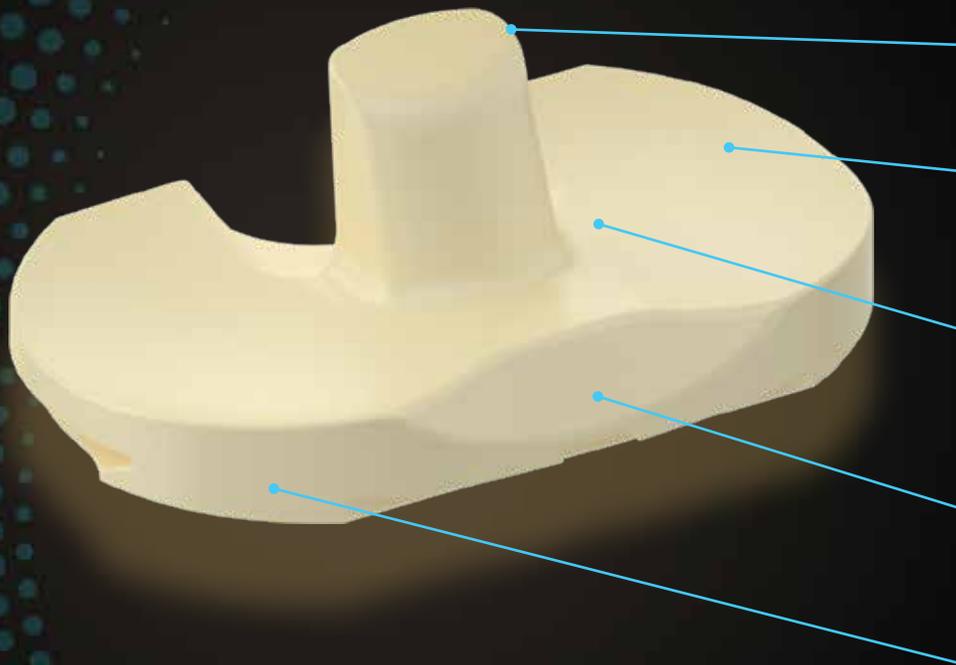

PS/PS+ Jump  
Height = 17mm

PS Constraint  
Free to Rotate

PS+ Constraint  
Constrained 3°  
IE and VV

Anterior Recess  
For Patella Ligament  
Wrapping


E-Dur Vitamin E  
Infused Polyethylene<sup>1</sup>

5° Sloped Post  
Increases Stability

Posterior Chamfer  
Designed to Reduce  
Impingement

0° Posterior Slope

Posterior and  
Central Dovetail  
Locking Feature

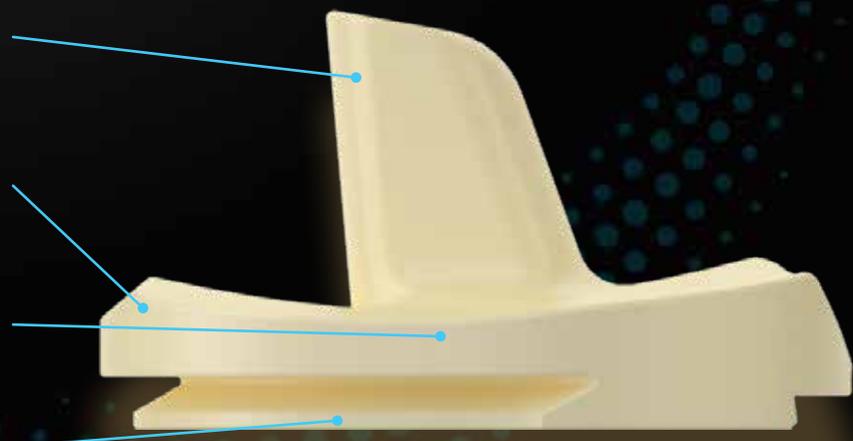

<sup>1</sup>Crosslinking w/ a 75 kGy gamma radiation dose in air then EtO sterilized.